


AirMini[™]



User guide English

$R_{\!\!X\,\text{Only}}$

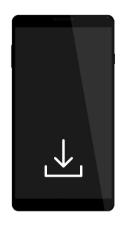









FNGLISH

Welcome

The AirMini™ system combines ResMed's AirMini self-adjusting pressure machine, masks and the AirMini app.



$oldsymbol{\Lambda}$ WARNING

Read this entire guide before using the machine.



In the US, Federal law restricts this device to sale by or on the order of a physician.

Indications for use

The AirMini self-adjusting system is indicated for the treatment of obstructive sleep apnea (OSA) in patients (female patients with mild to moderate OSA when using AutoSet for Her treatment mode) weighing more than 66 lb (30 kg).

It is intended for home and hospital use.

Contraindications

Positive airway pressure therapy may be contraindicated in some patients with the following pre-existing conditions:

- severe bullous lung disease
- pneumothorax
- pathologically low blood pressure
- dehydration
- cerebrospinal fluid leak, recent cranial surgery, or trauma.

Adverse effects

You should report unusual chest pain, severe headache, or increased breathlessness to your prescribing physician. An acute upper respiratory tract infection may require temporary discontinuation of treatment.

The following side effects may arise during the course of therapy with the machine:


- drying of the nose, mouth, or throat
- nosebleed
- bloating
- · ear or sinus discomfort
- eye irritation
- skin rashes.

At a glance

The AirMini system includes the following:

- AirMini: the AirMini machine, power supply unit, drawstring bag and one of the following:
- AirMini setup pack F20: Air tubing, F20 connector OR
- AirMini setup pack N20: Air tubing, N20 connector, HumidX[™], HumidX Plus OR
- AirMini mask pack for P10: Air tubing, AirFit P10 for AirMini mask (M), P10 cushion (S), P10 cushion (L), HumidX, HumidX Plus, AirFit P10 headgear clips.

Your AirMini system

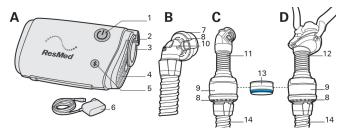

Your AirMini system is comprised of the AirMini machine and one mask pack or one setup pack.

Α	Airl	/lını	machine

- B F20 connector for AirMini and air tubing
- C N20 connector for AirMini and air tubing
- D AirFit P10 for AirMini and tubing
- 1 Start/Stop button
- 2 Power inlet
- 3 Air outlet
- 4 Air filter cover

Bluetooth[™] button

- 6 20W power supply unit
- 7 F20 connector for AirMini
- 8 Vent
- 9 Vent module
- 10 Anti Asphyxia Valve (AAV)
- 11 N20 connector for AirMini
- 12 AirFit P10 for AirMini
- 13 Heat moisture exchanger HumidX (blue) and HumidX Plus (gray)
- 14 AirMini tubing

Your machine is compatible with ResMed's F20 and N20 series masks and with the AirFit P10 for AirMini. Make sure all parts and accessories used with the machine are compatible. For compatibility information, refer to www.resmed.com for more information.

Your AirMini machine is designed to work with ResMed's AirMini app. However, the AirMini app is not essential to operate the machine. For more information on the AirMini app, refer to Using your machine with the AirMini app or the AirMini app Patient eHelp.

HumidX and HumidX Plus

The HumidX and HumidX Plus are waterless humidifiers. These are designed to provide more comfort by improving the moisture level in the air you breathe and offering relief from possible dryness.

The HumidX is a good starting point for finding the right level of humidification. The HumidX Plus will improve the humidity level further and is useful when the ambient humidity level is very low (eg, high altitudes, on an aircraft).

The HumidX and HumidX Plus can be used with the N20 connector for AirMini and AirFit P10 for AirMini only.

The HumidX and HumidX Plus must be replaced within 30 days of opening the HumidX/HumidX Plus pack.

Setting up and starting therapy

Using your mask

Before you set up your machine, make sure you use a compatible mask. For a full list of compatible masks for this device, see the Mask/Device Compatibility List on


www.resmed.com/downloads/masks. If you do not have internet access, please contact your ResMed representative.

For information on fitting your mask, refer to the mask user guide.

Use your AirMini machine only as directed by your physician or healthcare provider.

Set up your machine


- Plug the power supply unit into the machine and the power outlet. The connector has an angled edge so that the power cord can only be inserted into the machine one way.
- Connect the gray end of the AirMini tubing firmly to the air outlet.

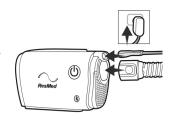

Add optional humidification (for N20 or AirFit P10 for AirMini only)


To insert the HumidX or HumidX Plus:

- 1. Open the connector by gently twisting it.
- Hold the HumidX/HumidX Plus at the sides, with the colored side pointing down and insert it
- 3. Gently push together the connector and twist until it locks in place.

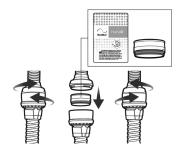

Connect your mask

- 1. Fit your mask as instructed in the mask user guide.
- 2. Attach the connector to your mask.

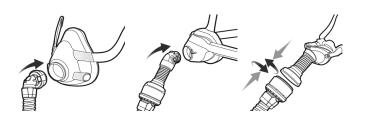

F20 Series

N20 Series

AirFit P10 for AirMini

Start therapy

 Start therapy by pressing the Start/Stop button or breathe normally if SmartStartTM is enabled.



Stop therapy

- 1. Remove your mask.
- Press the Start/Stop button or wait until the machine stops if SmartStop is activated.

Comfort features

Ramp Time, Pressure Relief and SmartStart are enabled on your machine by default. For more information or to change these settings, refer to the AirMini app Patient eHelp.

Using your machine with the AirMini app

Your AirMini machine can be used with ResMed's AirMini app. Once you have paired and connected your machine to your smart device via Bluetooth, you can start and stop therapy, change comfort settings and view your therapy data via the AirMini app. Data sent via Bluetooth is encrypted to protect your privacy.

Connecting your machine and smart device

Before you start, ensure that the latest version of the AirMini app is installed on your smart device. If not, download the app from the App Store $^{\text{TM}}$ or Google Play $^{\text{TM}}$.

- Ensure your machine is set up correctly and plugged into a power source.
- 2. On your smart device, enable Bluetooth.
- 3. Open the AirMini app.

The first time you open the AirMini app you will be asked to consent to the use of cookies and to accept ResMed's Terms and Conditions and Privacy Policy.

- 4. Check the boxes on the Welcome screen and then tap Continue.
- 5. On your machine, press the Bluetooth button 3 when prompted by the AirMini app.

When the light starts to flash blue \$\ \bar{\pi}\$, Bluetooth is enabled and your machine is ready to connect.

- 6. On your smart device, tap Connect.
 - When your machine is discovered, the machine name appears in the machine selection list.
- 7. Select the machine name from the list to connect.
- 8. The first time you pair your machine with your smart device, you will need to perform the authentication procedure.
 - When prompted by the AirMini app, authenticate by entering the four-digit key located on the back of your machine.

Alternatively, you can tap **Scan code**. You will be prompted to allow the use of your smart device's camera and then to position the machine's QR code within the camera's viewing frame. The QR code is found on the back of your machine.





9. Tap Continue to close the window.

When successfully connected, the Bluetooth connection icon is displayed in the top right corner of the app.

Your machine is in discoverable mode and

Bluetooth status indicators

Rhietooth indicator light is flashing.

*	blue.	ready for pairing.
*	Bluetooth indicator light is constantly blue.	Your machine is connected to your smartphone.
° *	Bluetooth indicator light is constantly white (airplane mode).	Bluetooth on your machine is disabled.
	Bluetooth indicator light is off.	Your machine is not connected to your smartphone.

Note: Once therapy starts, the indicator light will start to fade.

Cleaning and maintenance

It is important that you regularly clean your AirMini machine to make sure you receive optimal therapy. The following sections will help you with disassembly, cleaning and checking your machine.

⚠ WARNING

- As part of good hygiene, always follow cleaning instructions.
 Some cleaning products may damage the air tubing and affect its function, or leave harmful residual vapours that could be inhaled if not rinsed thoroughly.
- The AirMini connectors and AirFit P10 for AirMini contain a vent system to prevent carbon dioxide build-up in the mask. The vent system must be kept clean and clear of contamination in order to function correctly. Blocking or modifying the vent system may result in excessive carbon dioxide rebreathing.
- Regularly clean your AirMini tubing, connectors and mask components to receive optimal therapy and to prevent the growth of germs that can adversely affect your health.
- Regularly inspect the HumidX/HumidX Plus and follow the cleaning and maintenance instructions to prevent the growth of germs that can adversely affect your health.


△ CAUTION

If any visible deterioration of a system component is apparent (cracking, discoloration, tears etc.), the component should be discarded and replaced.

Disassembling

Remove the connector (F20 and N20 only)

 Remove the connector from the mask by pressing the side buttons and pulling it away

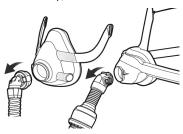

Disconnect the AirMini tubing (F20 only)

Disconnect the F20 connector from the AirMini tubing by gently twisting it and pulling it apart



Disconnect the AirMini tubing (N20 and AirFitP10 for AirMini only)


 Disconnect the vent module from the AirMini tubing by gently twisting it and pulling it apart.



Remove the HumidX/HumidX Plus (N20 and AirFit P10 for AirMini only)


 Disconnect the mask tube from the vent module by gently twisting it and pulling apart. If you are using HumidX or HumidX Plus, remove it from the vent module.

Note: The HumidX/HumidX Plus cannot be washed.

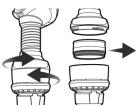

Unplug the machine

- Disconnect the AirMini tubing from the machine by pressing the side buttons on the cuff and pulling it away.
- Unplug the power supply unit from the power outlet and the machine.

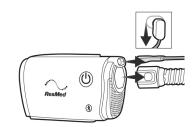

For mask disassembly instructions, see your mask user guide.

Cleaning and replacing parts

Detach all components according to the disassembly instructions and make sure that the HumidX/HumidX Plus is removed and stored in a clean, dry location.

Inspect the HumidX/HumidX Plus

- Check daily for any signs of damage or blockages caused by dirt or dust
- 2. It must be replaced within 30 days after opening.
- 3. When the HumidX/HumidX Plus is not being used, store it in a clean, dry location.

Note: HumidX/HumidX Plus cannot be washed.

Clean the AirMini tubing


- Clean and check the air tubing each week. AirMini tubing should be replaced at least every six months.
- 2. Wash the AirMini tubing in warm water using a mild liquid detergent, for example Ultra Joy dishwashing liquid.
- 3. Rinse thoroughly and allow to dry out of direct sunlight and/or heat.
- 4. Check the AirMini tubing and replace it if there are any holes, tears or cracks.

Note: Do not wash the AirMini tubing in a dishwasher or washing machine.

Clean the F20 connector for AirMini

- Clean and check the F20 connector daily. The connector should be replaced at least every six months.
- Wash in warm water using a mild liquid detergent, for example Ultra Joy dishwashing liquid.
- 3. Clean with a soft bristle brush paying particular attention to the vent holes.
- 4. Rinse well under running water. Shake to remove excess water and allow to dry out of direct sunlight.
- 5. Check that there is no dirt or dust in the vent holes.


Note: Do not wash in a dishwasher or washing machine.

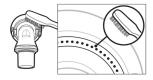

Clean the N20 connector for AirMini

- Clean and check the N20 connector daily. The connector should be replaced at least every six months.
- 2. Wash the connector in warm water using a mild liquid detergent, for example Ultra Joy dishwashing liquid.
- 3. Clean with a soft bristle brush paying particular attention to the vent holes inside the connector.
- 4. Rinse well under running water. Shake to remove excess water and allow to dry out of direct sunlight.
- 5. Check that there is no dirt or dust in the vent holes.

Note: Do not wash in a dishwasher or washing machine.

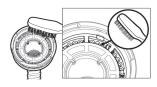

Clean the AirFit P10 for AirMini

Refer to your mask guide for cleaning instructions.


Change the air filter

Replace at least every 6 months. It cannot be washed.

Check and replace it more often if there are any holes or blockages caused by dirt or dust.

To change the air filter:

- 1. Slide the air filter cover off and remove the old air filter.
- 2. Insert a new air filter and replace the air filter cover.

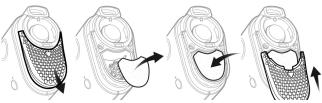

Note:

- Make sure the air filter cover is fitted at all times to prevent water and dust from entering the machine.
- The use of ResMed approved hypoallergenic filter will result in a small reduction in the accuracy of the delivered pressure at high leaks.

Cleaning the machine

1. Wipe the exterior with a dry cloth each week.

Reprocessing

The AirMini tubing, AirMini connectors, vent module and HumidX/HumidX Plus must be discarded and replaced, as they can not be disinfected and can not be used between other patients.

For mask cleaning instructions, see your mask user guide.

Traveling

You can take your AirMini machine with you wherever you go. Make sure you have the appropriate power supply for the region you are traveling to. For information on purchasing, contact your care provider.

Traveling by airplane

For some airlines, medical devices do not count toward carry-on luggage limits. Please check with your airline for their policy regarding medical equipment.

You can use your AirMini machine on a plane as it meets the Federal Aviation Administration (FAA) requirements. Air travel compliance letters can be downloaded and printed from www.resmed.com.

When using the machine on an airplane:

- When connected to power, disable Bluetooth (enter airplane mode) by pressing and holding the Bluetooth button ★ for at least ten seconds until the indicator light is white.
- Do not use the AirMini app.
- Use the Start/Stop button on your machine to start therapy,
- To reconnect Bluetooth (exit airplane mode), press the Bluetooth button.

Troubleshooting

If you have any problems, have a look at the following troubleshooting topics. If you are not able to fix the problem, contact your care provider or ResMed. Do not try to open the machine.

Problem/possible cause	Solution	
Air is leaking from around my mask / Mask is too noisy		
Mask may be fitted incorrectly.	Make sure your mask is fitted correctly. See your mask user guide for fitting instructions or run the Mask Fit function using the AirMini app to check your mask fit and seal.	

Problem/possible cause	Solution
The AirMini system may be assembled incorrectly.	Make sure your AirMini system is correctly assembled. Refer to Setting up and starting therapy for more information.
Air pressure in my mask seems too low	/ It feels like I am not getting enough air
Ramp may be in progress.	Using the AirMini app, confirm the Ramp function has been enabled. If so, wait for air pressure to build up or turn Ramp Time off using the AirMini app.
HumidX or vent might be blocked.	Inspect the HumidX and vent for blockages or damage. Refer to Cleaning and maintenance for more information.
HumidX may be wet.	The HumidX should be dry when starting therapy. Check if HumidX is wet and replace as necessary.
I cannot start therapy	
Power may not be connected.	Connect the power supply unit and make sure the plug is fully inserted. The green indicator light on the power supply unit should be illuminated.
	Check that the green indicator light above the Start/Stop button on the machine is also illuminated.
My therapy has stopped	
The air tubing may be disconnected.	Ensure that the air tubing is properly connected. Press the Start/Stop button to continue therapy.
You may have high leak and SmartStop	Make sure your mask is fitted correctly.
enabled.	Alternatively, disable SmartStop.

Problem/possible cause	Solution	
The Bluetooth indicator light is constantly white. I cannot pair to my machine		
Bluetooth is disabled.	Press the Bluetooth button for 3 seconds. When the Bluetooth indicator light flashes blue, the machine is ready for pairing.	
I cannot scan the QR code to pair my n	nachine	
The camera is not focused or broken, or the QR code label is damaged.	You can pair your machine manually by entering into your smart device the four- digit key found on the back of your machine.	
The QR code label is unreadable.	Contact your care provider.	
The light above the Start/Stop button is flashing green		
An error has occurred on the machine.	Unplug your machine. Wait a few seconds and then plug it back in.	
	If the problem persists, contact your care provider. Do not open the machine.	
I am experiencing dryness (dry or blocked nose)		
Mask may be fitted incorrectly	Adjust and refit your mask to improve seal. Check that there is no air leaking from the mask.	
	If mask seal is good, try the HumidX Plus.	
Humidification may be required	Try therapy using the HumidX	

Problem/possible cause	Solution	
I am using the HumidX/HumidX Plus and I am getting uncomfortable droplets of water on my nose, mask or air tubing		
Humidity level is too high	If using the HumidX: Remove the HumidX and try therapy without humidification.	
	If using the HumidX Plus: Try the Humid X instead. This will help to reduce the water droplets.	
	If you are still experiencing problems, contact your care provider.	
	Note: You may need to use the HumidX/HumidX Plus when ambient humidity conditions change.	

General warnings and cautions



- Make sure that you arrange the air tubing so that it will not twist around the head or neck.
- · Keep the power cord away from hot surfaces.
- Make sure the power cord and plug are in good condition and the equipment is not damaged.
- If you notice any unexplained changes in the performance of the machine, if it is making unusual sounds, if the machine or the power supply are dropped or mishandled, or if the enclosure is broken, discontinue use and contact your care provider or your ResMed Service Center.
- Do not open or modify the machine. There are no user serviceable parts inside. Repairs and servicing should only be performed by an authorized ResMed service agent.

Δ warning

- Beware of electrocution. Do not immerse the machine, power supply or power cord in water. If liquids are spilled into or onto the machine, unplug the machine and let the parts dry. Always unplug the machine before cleaning and make sure that all parts are dry before plugging it back in.
- Do not use with oxygen. Any sources of oxygen must be located more than 1 m away from the machine to avoid risk of fire and burns.
- Do not perform any maintenance tasks while the machine is in operation.
- Use only AirMini masks and accessories that allow for normal breathing to avoid asphyxiation.
- The use of accessories other than those specified for the machine is not recommended. This may result in increased electromagnetic emissions or decreased immunity of the machine and result in improper operation.
- The machine should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, the machine should be observed to verify normal operation in the configuration in which it will be used.
- The F20, N20 and P10 connectors for AirMini are fitted with vents. The F20 connector for AirMini is also fitted with an AAV. The AAV and vents have specific safety functions to prevent carbon dioxide build-up in the mask. The AAV and vents should be kept clean and not covered or blocked. The connectors should not be used if the vent or AAV are damaged as they will not be able to perform their safety functions. The connectors should be replaced if the valves of the vent or AAV are damaged, distorted or torn.
- The machine is not intended to be operated by persons (including children) with reduced physical, sensory or mental capabilities without adequate supervision by a person responsible for the patient's safety.

 The device has not been tested or certified for use in the vicinity of X-ray, CT or MRI equipment. Do not bring the device within 13 ft (4 m) of X-ray or CT equipment. Never bring the device into an MR environment.

⚠ CAUTION

- Use only ResMed AirMini parts, masks and accessories with the machine. Non-ResMed AirMini parts may reduce the effectiveness of the treatment, may result in excess carbon dioxide rebreathing and/or damage the machine. For compatibility information, refer to www.resmed.com for more information.
- Blocking the air tubing and/or air inlet of the machine while in operation could lead to overheating of the machine.
- Keep the area around the machine dry, clean and clear of anything (eg, clothes, pillows or bedding) that could block the air inlet or cover the power supply unit.
- Do not use bleach, chlorine, alcohol, or aromatic-based solutions, moisturizing or antibacterial soaps or scented oils to clean the machine or air tubing. These solutions may cause damage and reduce the life of the products.
- Do not insert any USB cable into the machine or attempt to plug the power supply unit into a USB device. This may cause damage to the machine or USB device.

Technical specifications

Units are expressed in cm H₂O and hPa. 1 cm H₂O is equal to 0.98 hPa.

Power supply unit

AC input range: 100–240V, 50–60Hz

115V, 400Hz (nominal for aircraft use)

DC output: 24V ______ 0.83A

Typical power consumption: 6.3W Peak power consumption: 27W

If power is interrupted during therapy, the machine will automatically restart therapy when power is restored.

AirMini power	
Standby power consumption:	1W
Environmental conditions	
Operating temperature:	+41°F to +95°F (+5°C to +35°C)
	Note: Under extreme ambient temperature conditions (104°F/40°C) the air and parts of the tubing near the mask may reach temperatures up to 109°F/43°C. These parts of the tubing are not in contact with the user and the machine remains safe under these extreme conditions.
Operating humidity:	10 to 95% relative humidity, non-condensing
Operating altitude:	Sea level to 8,500' (2,591 m); air pressure range 1013 hPa to 738 hPa
Storage and transport temperature:	-13°F to + 158°F (-25°C to + 70°C)
Storage and transport humidity:	5 to 95% relative humidity, non-condensing

Electromagnetic compatibility

The AirMini complies with all applicable electromagnetic compatibility requirements (EMC) according to IEC60601-1-2:2014, for residential, commercial and light industry environments. Portable and mobile RF communications equipment should be used no closer to any part of the machine, including cables, than the recommended 10 cm separation distance. The AirMini has been designed to meet EMC standards. However, should you suspect that the machine performance (eg, pressure or flow) is affected by other equipment, move the machine away from the possible cause of interference.

The AirMini complies with Part 15 of the FCC Rules and Industry Canada license-exempt RSS standards. Operation is subject to the following two conditions: This machine may not cause harmful interference, and this machine must accept any interference received, including interference that may cause undesired operation.

FCC ID: QOQBT121, IC: 5123A-BGTBT121

Additional information regarding the FCC Rules and IC compliance for this machine can be found on www.resmed.com/downloads/devices.

IEC 60601-1 (Edition 3.1) classification

Class II (double insulation), Type BF, Ingress protection IP22.

Maximum single fault steady pressure

Machine will shut down in the presence of a single fault if the steady state pressure exceeds 30 cm H_2O (30 hPa) for more than 6 sec or 40 cm H_2O (40 hPa) for more than 1 sec.

Sound

Declared dual-number noise emission values in accordance with ISO 4871:1996

Pressure level measured according to ISO 30 dBA with uncertainty of 2 dBA

80601-2-70:2015 (CPAP mode)

Power level measured according to ISO 80601- 38 dBA with uncertainty of 2 dBA

2-70:2015 (CPAP mode)

The A-weighted sound pressure level measured according to ISO 17510:2015:

Air tubing connected to AirMini F20 connector 19 dBA with uncertainty of 3 dBA with mask

Air tubing connected to AirMini N20 or P10 17 dBA with uncertainty of 3 dBA

connectors with mask

The A-weighted sound power level measured according to ISO 17510:2015:

Air tubing connected to AirMini F20 connector 27 dBA with uncertainty of 3 dBA

with mask

Air tubing connected to the AirMini N20 or P10 25 dBA with uncertainty of 3 dBA

connectors with mask

AirMini machine

Dimensions: 5.4" (W) x 3.3" (D) x 2.1"(H)

136 mm (W) x 84 mm (D) x 52 mm (H)

Weight: 300g (10.6oz)

Housing construction: Flame retardant engineering thermoplastic

Proprietry connector, 16mm inner

diameter. Not compatible with EN ISO 5356-1 connectors

Air filter

Air outlet:

Standard: Material: Polyester non woven fiber

Average arrestance: >75% for ~7 micron

dust

Hypoallergenic: Material: Acrylic and polypropylene fibers

in a polypropylene carrier

Efficiency: >98% for ~7-8 micron dust;

>80% for ~0.5 micron dust

AirMini air tubing

Material: Flexible plastic
Length: 6'2" (1.9 m)
Inner diameter: 0.6" (15 mm)

Aircraft use

ResMed confirms that the machine meets the Federal Aviation Administration (FAA) requirements (RTCA/D0-160, section 21, category M) for all phases of air travel.

Wireless technology

Technology used: Bluetooth

Connection types: SPP, iAP2, GATT Frequency: 2402 to 2480 MHz

 $\begin{array}{ll} \text{Max RF power output:} & +4 \text{ dBm} \\ \text{Operating range:} & 10 \text{ m (Class 2)} \end{array}$

It is recommended that the machine is at a minimum distance of 0.43" (1.1cm) from the body during operation. Not applicable to masks, air tubing or accessories.

Declaration of Conformity to 1999/5/EC

(DoC to the R&TTE Directive)

ResMed declares that the AirMini device is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. A copy of the declaration of conformity (DoC) can be found on www.resmed.com/downloads/devices.


Operating pressure range

AutoSet, AutoSet For Her, CPAP: 4 to 20 cm H₂O (4 to 20 hPa)

Design life

Machine, power supply unit: 5 years
Air tubing, connectors: 6 months
HumidX: 30 days

Pneumatic flow path

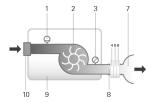


- 1. Flow sensor
- 2. Blower
- 3. Pressure sensor
- 4. AAV (F20 connector only)
- 5. Vent
- 6. HumidX (N20, P10 connectors only)
- 7. Mask
- 8. Air tubing
- 9. Machine
- 10. Inlet filter

General

The patient is an intended operator.

Pressure flow curve (ISO17510:2015)

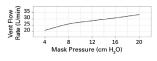

Pressure [cm H ₂ O (hPa)]	Flow [L/min]
4	21
8	26
12	29
16	31
20	33

Inspiratory and expiratory resistance with AAV open to atmosphere (ISO17510:2015)	F20 connector
Inspiration at 50L/min	0.7 cm H₂O (hPa)
Expiration at 50L/min	0.5 cm H₂O (hPa)
AAV Pressures (ISO17510:2015)	F20 connector
Open-to-atmosphere pressure	0.9 cm H₂O (hPa)
Closed-to-atmosphere pressure	1.0 cm H₂O (hPa)

Measurement system uncertainties

In accordance with ISO 80601-2-70:2015 the measurement uncertainty of the manufacturer's test equipment is:

For measures of flow	\pm 1.5 L/min or \pm 2.7% of reading (whichever is greater)
For measures of static pressure	$\pm 0.15 \text{ cm H}_2\text{O (hPA)}$
For measures of dynamic pressure	$\pm 0.27 \text{ cm H}_2\text{O (hPA)}$
For measures of volume (< 100 mL)	± 5 mL or 6% of reading (whichever is greater)
For measures of volume (≥ 100 mL)	± 20 mL or 3% of reading (whichever is greater)
For measures of time	± 10 ms

Note: ISO 80601-2-70:2015 stated accuracies and test results provided in this manual for these items already include the relevant measurement uncertainty from the table above.

Pressure accuracy

Maximum static pressure variation at 10 cm H_2O (10 hPa) according to ISO 80601-2-70:2015:

Tested with full face mask: ±0.5 cm H₂0 (0.5 hPa)

Maximum dynamic pressure variation according to ISO 80601-2-70:2015

Machine with full face mask

Pressure [cm H₂O (hPa)]	10 BPM	15 BPM	20 BPM
4	0.5	0.7	1.0
8	0.5	0.7	1.0
12	0.5	0.7	1.0
16	0.5	0.7	1.0
20	0.5	0.7	1.0

Flow (maximum) at set pressures

The following are measured accordingly to ISO 80601-2-70:2015 at the end of the specified air tubing:


Pressure cm H ₂ O (hPa)	AirMini machine and AirMini air tubing L/min
4	119
8	116
12	112
16	108
20	105

Notes:

 The manufacturer reserves the right to change these specifications without notice.

The mask system does not contain PVC, DEHP or phthalates.This product is not made with natural rubber latex.
Symbols
The following symbols may appear on the product or packaging.
Read instructions before use. Indicates a warning or caution.
Follow instructions before use. Manufacturer. EC REP European Authorized Representative. EC REP Batch code REF Catalog number. SN Serial number. DN Device number. On / Off. IP22 Protected against finger sized objects and against
dripping water when tilted up to 15 degrees from specified orientation. Direct current. Type BF applied part. Class II equipment. Humidity limitation. Temperature limitation. Non-ionising radiation. Bluetooth. Rx Only Prescription only (In the US, Federal law restricts these devices to sale by or on the order

of a physician). Operating altitude. Atmospheric pressure limitation. Complies with RTCA DO-160 section 21, category M.

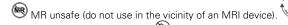

Replace 30 days after opening. FRACKAGE Do not use if package is damaged.



Environmental information

This machine should be disposed of separately, not as unsorted municipal waste. To dispose of your machine, you should use appropriate collection, reuse and recycling systems available in your region. The use of these collection, reuse and recycling systems is designed to reduce pressure on natural resources and prevent hazardous substances from damaging the environment.

If you need information on these disposal systems, please contact your local waste administration. The crossed-bin symbol invites you to use these disposal systems. If you require information on collection and disposal of your ResMed machine please contact your ResMed office, local distributor or go to www.resmed.com/environment.

Servicing

The AirMini machine is intended to provide safe and reliable operation when operated in accordance with the instructions provided by ResMed. ResMed recommends that the AirMini machine be inspected and serviced by an authorized ResMed Service Centre if there is any sign of wear or concern with device function. Otherwise, service and inspection of the products generally should not be required during their design life.

Limited warranty

ResMed Ltd (hereafter 'ResMed') warrants that your ResMed product shall be free from defects in material and workmanship from the date of purchase for the period specified below.

Product		vvarranty period
•	Mask systems (including mask frame, cushion, headgear and tubing)—excluding single-use devices	90 days
•	Accessories—excluding single-use devices	
•	Batteries for use in ResMed internal and external battery systems	6 months
•	CPAP (including external power supply units)	2 years

Marrantunariad

This warranty is only available to the initial consumer. It is not transferable.

If the product fails under conditions of normal use, ResMed will repair or replace, at its option, the defective product or any of its components.

This Limited Warranty does not cover: a) any damage caused as a result of improper use, abuse, modification or alteration of the product; b) repairs carried out by any service organization that has not been expressly authorized by ResMed to perform such repairs; and c) any damage or contamination due to cigarette, pipe, cigar or other smoke.

Warranty is void on product sold, or resold, outside the region of original purchase.

Warranty claims on defective product must be made by the initial consumer at the point of purchase.

This warranty replaces all other expressed or implied warranties, including any implied warranty of merchantability or fitness for a particular purpose. Some regions or states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

ResMed shall not be responsible for any incidental or consequential damages claimed to have resulted from the sale, installation or use of any ResMed product. Some regions or states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

Dradust

This warranty gives you specific legal rights, and you may also have other rights which vary from region to region. For further information on your warranty rights, contact your local ResMed dealer or ResMed office.

Further information

If you have any questions or require additional information on how to use the machine, contact your care provider.









ResMed Ltd I

Distributed by

1 Elizabeth Macarthur ResMed Corp 9001 Spectrum Center Boulevard San Diego CA 92123 Drive Bella Vista NSW USA [ECIREP] ResMed (UK) Ltd 96 Jubilee Ave Milton Park Abingdon 2153 Australia Oxfordshire OX14 4RW UK

See ResMed.com for other ResMed locations worldwide. AirMini, AutoSet, HumidX and SmartStart are trademarks and/or registered trademarks of the ResMed family of companies. For patent and other intellectual property information, see ResMed.com/jp. "Made for iPhone" means that an electronic accessory has been designed to connect specifically to iPhone and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance. Apple, the Apple logo, iPhone and App Store are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. IOS is a trademark of Cisco, licensed to Apple Inc. Bluetooth is a registered trademark of the Bluetooth SIG, Inc. Joy is a trademark of Procter & Gamble. The Android robot is reproduced or modified from work created and shared by Google and used according to terms described in the Creative Commons 3.0 Attribution License. Android, Google Play and the Google Play logo are trademarks of Google Inc. Galaxy, Edge and Note are trademarks of Samsung. This product uses software developed by the ANTIR Broiet (http://www.antir.) provided in 388010/1 2017-01

200010